



BENEFITS & FEATURES



TIME-EFFICIENT

Hearing screening in less than 2 minutes.



CUSTOMISABLE PASS/REFER CRITERIA

Default and customisable settings based on screening frequencies and intensities.



ADJUSTABLE PROTOCOLS

Customisable screening intensities for adults and children.



THRESHOLD TRACKING

Optional re-screen and severity screening for threshold tracking to determine hearing loss classification.



ENVIRONMENTAL NOISE MONITORING

Real-time monitoring of ambient noise to alert the user of environmental noise concerns.



COST-EFFECTIVE

Decrease screening costs by >50% and perform more tests in less time.



CLINICALLY VALID TESTS

Evidence-based with calibrated headphones.



DIGITAL DATA MANAGEMENT

Patient, test and facility data consolidated instantly on a secure online database, mHealth Studio Cloud.

"hearScreen® is an excellent low-cost smartphone-based screening audiometer. It enabled our audiology department to provide more efficient hearing screening services to school children and adults." - **MESHACK MOEPENG**, Audiologist, Bamalete Lutheran Hospital


TECHNICAL SPECIFICATIONS

Pure tone test frequencies | 0.5 to 8 kHz

Rise/fall time | ~35 ms

Total harmonic distortion | <2%

Test intensity levels | 20 to 70 dB HL

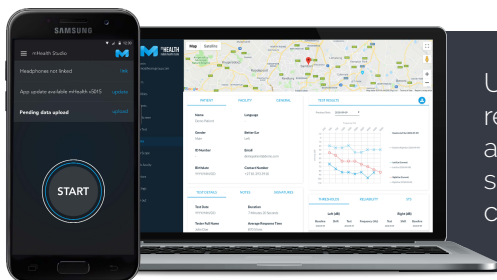

Use this secure cloud portal to enable remote support, patient counselling, assistance and reporting. All test data seamlessly integrates with the secure cloud-based data management system.

MHEALTH STUDIO CLOUD

FDA Registered | CE Certified | TGA Certified | ISO 13484 Certified | SAHPRA Registered | HIPAA Self-declared | POPIA Self-declared | GDPR Self-declared

+1 415 212 5500 | +27 (0)12 030 0268 | info@hearXgroup.com | www.hearXgroup.com



@hearXgroup



@hearX_Group



hearX Group



@hearXgroup

© 2019-2021 hearX IP (Pty) Ltd. Copyright in this work including its layout and the designs contained herein is reserved and the making of copies is prohibited.